

Astell Scientific Ltd

19-21 Powerscroft Road

Sidcup

Kent

DA14 5DT

United Kingdom

Tel: +44 (0)20 8300 4311

Fax: +44 (0)20 8300 2247

email: info@astell.com

Website: astell.com

Product Specification

Model Ref.

AMA272BT



Astell AMA272BT – 130 litre Top Loading Autoclave – DIRECT STEAM VERSION

Full Colour Touch Screen Controller

Ideal for Liquids & Most laboratory applications, including Discard, Glassware, Laboratory Instruments etc

Standard Features

All the following are Value - Added Features
Included AS STANDARD

- Full Colour Touch Screen Controller
- 5 Programs
- Swiftlock rapid access safety closure
- Over temperature water cut out
- Thermocouple entry port
- Low Water Indicator
- Automatic timed air purge system
- Cool Lock Protection (fluid Cycles)
- Pressure gauge
- Vent valve
- Electro Polished Stainless Steel Pressure Vessel
- Safety Valve Testing
- On Screen Graph display
- Simple connection to on-site steam



Options

Integral Data printer
10 Program Controller
Datalogger
Load Sensed Process Timing
Advanced Water Cooling

RS232
Ethernet
Assisted Air Cooling
AutoFill

Accessories

½ Depth Stainless Steel Basket
Full Depth Stainless Steel Basket
½ Depth Morrison Discard Container
Full Depth Morrison Discard Container
Load Support Plate

See options & Accessory pages for full details

Overall Capacity	Chamber Diameter	Working Depth	Overall Dimensions	Operating Range
130 Litres	456mm	680mm	660x840x1150 mm (wxdxh)	100 -138°C (0.2-2.4 bar)

Power Requirements	Single Phase 2A (Steam)
---------------------------	--------------------------------

Steam & Drainage Requirements	Steam: A supply of dry saturated steam at 4 bar max Drainage: Floor level, 35mm, ideally non-manifolded, to withstand freeflowing steam
------------------------------------------	--------------------------------------------------------------------------------------------------------------------------------------------

Approx. Nett Weight:	Shipping Weight	Shipping Dimensions
242 Kg	268 Kg	109x79x151cm (wxdxh)

Options & Accessories

The basic model quoted above is manufactured to the highest standards. It will be fully calibrated and tested prior to dispatch, and will be suitable for a wide range of sterilization procedures.

Please note however, that the following factory-fitted options are available for the above model. For certain applications, these options can provide increased flexibility and improved performance. If purchased with the autoclave, there would be no additional freight/FOB costs incurred.

AAP800 Multi Program Controller: A multi-program controller which allows storage of up to 10 programs for added speed and versatility in sterilization procedures. Additional features on this option include a digital pressure display and a "menu" allowing the user to display and cross check all programs at the press of a button. If the optional data printer (AAR100) is fitted, this may be printed out.

Data Monitoring & Recording Options

AAR100 Integral Data Printer: This printer provides a permanent and traceable record of: "Time "Temperature "Pressure "Batch No. "Cycle Name. There is also provision for operator signature; the printer also provides reports of cycle settings and servicing information.

AAR120 RS232 Interface: This allows cycle progress to be monitored on an external computer equipped with suitable software.

AAR122 Ethernet Interface: Via VNC Viewer. Allows visual display of autoclave touch screen

Important option to be considered if fluids are being processed

AAN014 Load Sensed Process Timing: Load Sensed Process Timing allows the sterilization cycle to be controlled via the temperature achieved in the centre of the load. The operation of load sensed process timing is controlled via the 'Secure' programmer. A 'wandering' probe, situated within the chamber is inserted into the load, or load simulator, and initiates the sterilization period once the probe reaches the programmed threshold temperature. A selectable 'Profiled Overshoot Boost' speeds up the cycle and minimises over-processing of media loads.

Cooling Options

AAP006 Assisted Air Cooling: (Non Vacuum Models) A powerful fan-assisted cooling system which greatly reduces cool-down time. The operation of the Assisted Air Cooling system is controlled via the 'Secure' programmer and is available on all models.

AAP080 Advanced water cooling Water is circulated through cooling coils in contact with chamber, resulting in rapid decrease of temperature – advanced water cooling offers substantial time saving and is ideal for the busy Laboratory.

Additional Options

AVC100 Exhaust Heat Exchange System: An option to significantly reduce the temperature of the exhaust/waste in situations where heat resistant drains are not present. N.B. This option requires a mains water supply.

AAN009 Pulsar Free Steaming: (Non Vacuum Models) Selected via the 'Secure' Controller, this optional feature ensure adequate air removal from 'difficult' loads. The cyclical opening and closing of the vent valve greatly increases turbulence at the Freesteaming stage. Pulsar freesteaming All non-vacuum models

AVC005 Option For CAT III Compliance - Non Vacuum Models

For models without vacuum, or with heaters in the chamber, the bacterial retentive filter is fitted to the exhaust line and is sited to ensure that it is sterilized each cycle

IQ/OQ Documentation Pack

IQ Documentation

Details of calibration equipment * Order Acknowledgement * PED (Pressure Equipment Directive) Compliance * Declaration of Conformity * FAT (Factory Acceptance Test) * Drawing Schedule * ISO 9001:2000 Certification * Pressure vessel specification * Door safety checks

OQ Documentation


Chamber temperature distribution (per cycle) * Automatic control test (per cycle)

Please note: This is our standard IQ/OQ Documentation package. If other documents are required, please contact us with details of your specific requirements.


Chart Recorders

Astell are able to offer a comprehensive range of Chart Recorders. Please either contact us or visit our website (www.astell.com) for further information

Accessories

- AAN036** Stainless steel wire basket(400dx400hmm) Top loading models
 - AAN042** Stainless steel wire basket(400dx210hmm) Top loading models
 - AAN048** Container with perforated base (380dx400hmm) Top loading models
 - AAN056** Morrison discard container (390dx500h mm) Top loading models
 - AAN058** Morrison discard container (390dx355h mm) Top loading models
- 

Full Colour Touch Screen Controller

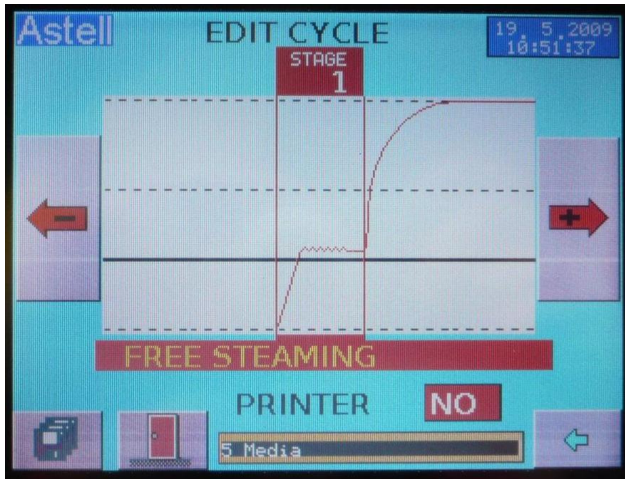

- Full Colour Touch Screen
- Data Archiving
- 5 Selectable security operating levels with Password Protected Security
- Easy Fingertip Operation even when wearing Lab Gloves
- Easy selection of program by Name and Settings
- Audible confirmation of key presses
- Simple Cycle Selection
- RS232 Interface / Port - Allowing connection to other peripheral devices such as a printer
- Operator Interface
- Built-in provision for mandatory Safety valve tests etc
- Enables continual monitoring

- Warnings of Faults
- Onscreen Display of Program Info, Times, Temperatures, Pressure, Cycle stage, Cycle Counter etc
- Internal Fault Detection with automatic fault diagnosis and complete fault history for easy servicing
- Simple to use Controls set out on the large illuminated Touch Screen display with menus for all operations
- Optional Communications including Ethernet and RS232

This control system is an advance in sterilization control technology, bringing together years of unrivalled experience, to produce a user friendly, fully automatic control system, to meet and exceed the expectations of even the most demanding laboratories and centres of sterilization.

The controller consists of a wipe clean touch screen based on an industrial PLC controller, combined with a number of analogue and digital input/output modules.

The controller software has been developed by Astell for the precision control of autoclaves.

The code is split into well-defined blocks, each of which performs a specific function and has been written in high level language to allow maximum clarity.

NB – This unit does not have a vacuum feature and is therefore unsuitable for wrapped instruments etc

Safety Features:

Design to PD97/23/EC

Overheat protection

Insulated safety door

Cooling locks

(in accordance with H.S.E. PM73) preventing opening of autoclave above 80°C. (for fluid & discard cycles)

Alarms

for Cycle Fault - Cycle Interruption - Sterilize Failure - Water Low - Door Unlocked

Door Seal

Self Energising / Service independent

Door

The door release is interlocked by both temperature and pressure to ensure all residual pressure has completely and effectively vented to atmosphere before the doors can be opened. The door will retain its positions in the event of failure of any service. The door is thermally insulated to prevent the surface temperature presenting a hazard to operators. The surface temperature will not exceed IEC 61010 requirements. A cycle cannot start until the door is closed and locked. Steam cannot be applied to the chamber unless the door is closed and locked.

Interlocks

Safety interlocks are provided, and are achieved by hardware, separate from and additional to the control system. All interlocks are configured to fail-safe and to provide a signal to the control system to indicate that normal operation has been prevented, and to terminate the current cycle. The interlock system is designed so that its function can be tested during routine maintenance. Safety related interlocks are either hard wired or piped. The following safety interlocks are provided: ♦If the door is not closed, the steam supply to the chamber will be isolated ♦If the pressure in the chamber exceeds 0.15 bar the door will remain locked.

Performance Tests

All electrical equipment is Safety Tested in accordance with the Low Voltage Directive.

Astell shall perform the following standard Factory Acceptance Tests. The tests are included in the machine costs as per the quotation prior to despatch:

All Astell autoclaves are fully tested and calibrated before despatch in line with our ISO9001-2000 procedures.

Connect and check all supplies	Produce printout for each cycle tested (When printer option fitted)
Check software version	Check door interlocks
Power up controller and check door open/close operation	Carry out sterilize monitor timer test on each cycle
Install relevant cycle data as required by client	Run cycle and check for any leaks
Calibrate all temperatures and pressures	Check all safety valves
Check rotation of all pumps and motors (if applicable)	Run each cycle and check conforms with applicable standards
Check and document safety devices	Document in Astell procedures

Applicable Standards

PED EN/97/23/EC /

ISO9001-2000 / UKAS / IEC 61010

Medical Devices Directive 93/42/EEC

Medical Devices Quality Management System - BS EN ISO 13485:2003

Autoclave Safety

N.B. Please note that all Astell autoclaves are manufactured to the highest standards and in full compliance with the Pressure Equipment Directive – i.e. PD5500/PED/97/23EC. Whilst all units have the necessary safety features to minimise user risk, and help ensure long term reliability, it is recommended that:

- a) Routine safety checks are carried out in accordance with Astell manuals and in compliance with current pressure regulations and/or insurance requirements.
- b) Autoclaves are serviced regularly by Astell or Astell trained/recommended engineers. *(UK only: Please contact us for further information and costs on our range of Preventative Maintenance contracts).*